




PA-249

**Force Table, Model PA-249**, it is used to verify the law of composition and resolution of forces. This consists of machined table 40cm in diameter. It raised rim has a scale, graduated in 360 degree, mounted on vertical support rod based provided with levelling screws. Body under study is a ring in center of the table fastened by cords passing over three pulleys, attached with the table. Complete with three pulleys, three weight hangers (50gm), each hanger carries 9 x 20gm., 1 x 10gm., 2 x 5gm.

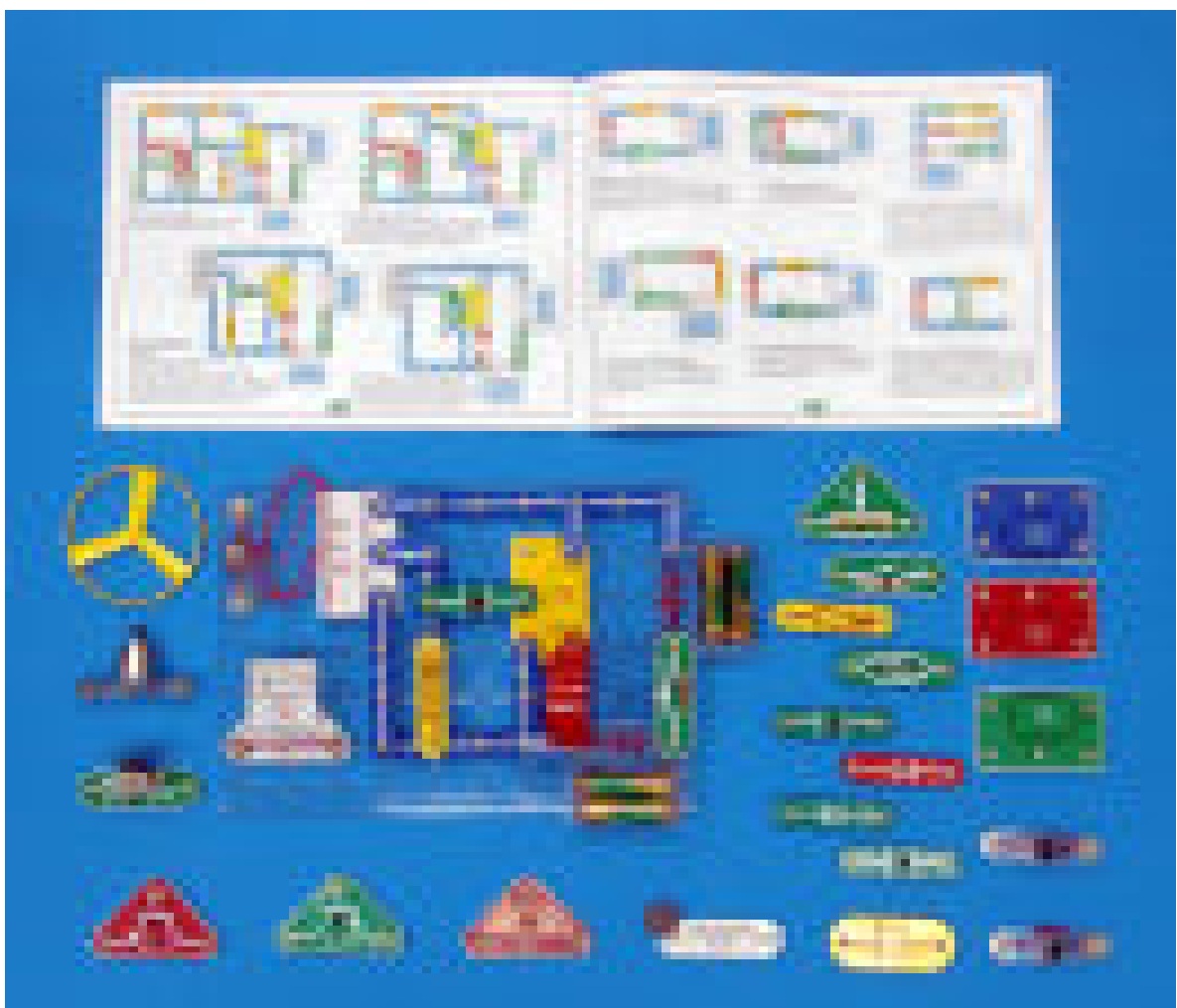

**Model PA-260**

**Model PA-260 Electronic Experiment Kit**, ideal for safely teaches the basic principles of electronics. The kit composed of a number of electronic component blocks and wires of different length, which can be used to build 688 different circuits, such as light-controlled buzzer, touch-controlled fan and more. Circuit block are made of plastic. With instruction.